

HerbaChaud®

The natural heat patch




www.herbachaud.com

Heat – Minerals – Herbs

How does the HerbaChaud® heat patch work?

The HerbaChaud® heat patch combines proven naturopathy methods with the millennia-old knowledge of traditional Chinese medicine.

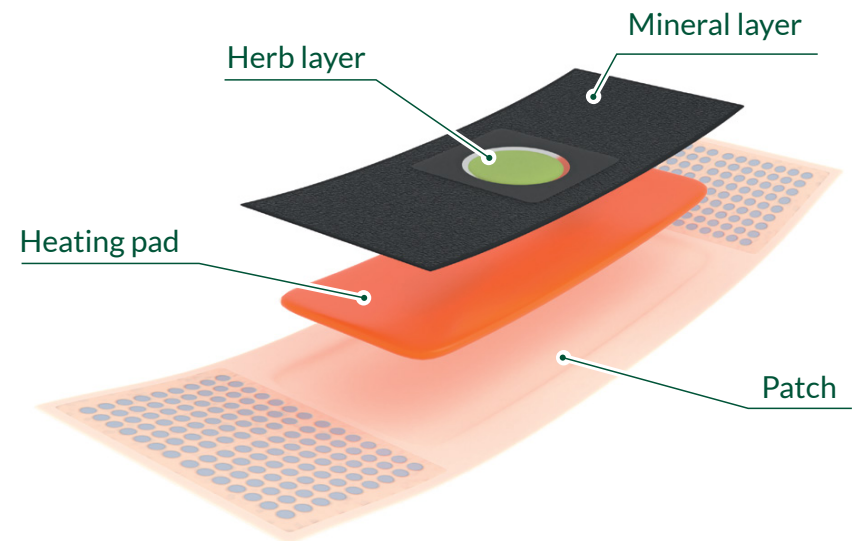

With the soothing effects of natural minerals, the healing power of traditional Chinese herbs and relaxing deep heat, HerbaChaud® heat patches bring together three active ingredients in a modern and easy-to-use application.

The heat patch produces a constant temperature of up to 40 °C directly at the pain trigger point, penetrating into the lowest muscle layers. This deep heat, lasting up to 12 hours, relaxes the muscles and specifically promotes blood circulation in the tissue. This improves the supply of oxygen and nutrients to the muscles, thereby promoting the natural self-healing process.

Stay active with the **HerbaChaud**[®] heat patch

For application on the back, neck, abdomen or knees in cases of:

- 🍃 backache and lower back pain
- 🍃 osteoarthritis and joint complaints
- 🍃 neck tension and shoulder area
- 🍃 strains and sprains
- 🍃 delayed onset muscle soreness (DOMS)
- 🍃 menstrual cramps



The **HerbaChaud**[®] heat patch is naturally activated by atmospheric oxygen. It is especially effective if applied at night: apply the **HerbaChaud**[®] heat patch to the painful area about 2 hours before going to bed. The intensive, deep heat unfolds its long-lasting, soothing and pain-relieving effect during the entire night's sleep.

HerbaChaud[®] is also suitable for mechanically stressed body parts, such as knees or elbows.

Trust the test winner!

Gesundheitstipp

Test Winner

Test Gesundheitstipp n° 2/2022

Over
93%

of users report
significant pain
relief.

HerbaChaud[®] has been selected as the test winner among thermal patches in several independent tests.

Try for yourself:

HerbaChaud[®] – the deep-acting heat patch that makes a noticeable and natural difference!



HerbaChaud[®]
The natural heat patch

The HerbaChaud[®]
heat patch helps
you to keep active,
naturally and
gently

The natural patch that heats the longest!

You will feel these benefits

- 🌿 100% natural plant and mineral active ingredients
- 🌿 completely free of side effects
- 🌿 no risk of skin irritations
- 🌿 also ideal for long-term use



suitable for people
with allergies



12 hours
of heat



100%
natural

The natural heat patch

HerbaChaud® natural heating patches are physically acting, odorless heating patches, which emit deep penetrating warmth (40°C) for up to 12 hours. The long-lasting warmth promotes blood circulation in muscles and joints.

APPLICATIONS

- Muscle and joint pain due to muscle tension, osteoarthritis, overexertion, strains, sprains
- Back pain/sacroiliac pain
- Tension in the neck and shoulders
- Menstrual cramps

HOW TO USE

The product is activated by contact with oxygen in the air. Only open immediately before use. **HerbaChaud®** is applied directly on the skin and can be worn under clothing. Apply the patch to the affected area.

First, make sure the skin is clean and dry. Open the foil packaging, then remove the two protective films and apply the patch to the affected area, adhesive side down. If your skin is sensitive, apply the patch instead with the adhesive facing the inside of your undergarment. The patch can also be fixed on the sides with medical adhesive tape. The patch will heat up after 30 minutes and stay warm for up to 12 hours. To ensure that enough oxygen passes through the back of the patch, do not obstruct the air supply with tight clothing. **This product is for single use only!**

INGREDIENTS

- Iron, coal, table salt, water, herbal mixture

RISKS

- For external use only. Do not swallow, contains substances harmful to your health. If swallowed, consult a doctor immediately.
- Do not heat in a microwave or try to heat the plaster a second time. The plaster could catch fire.

CONTRAINDICATIONS

- Do not use on face, chest, inner arms, armpits, fingers, toes, injured skin, bruises or swellings.
- Do not use on open wounds, near mucous membranes, the eyes or in the bend of the knee.
- Do not use in combination with other dermatics, heat sources or heat lotions/ ointments.

For people over 55 years of age: The risk of burns increases with age. When applying the heat patch to the lower back or hips, the patch should be worn over a layer of clothing and not directly on the skin. When applying over the knee, the heat patch should be worn over a thin cloth, such as a washcloth, and not applied directly to the skin. **Do not wear while sleeping.**

SAFETY INSTRUCTIONS

- For single use only
- For external use only
- Do not ingest, as the substances are harmful to health.

If accidentally swallowed, seek medical attention immediately.

- Do not heat in the microwave or attempt to heat the patch. There is a risk of the patch catching fire.
- All heat patches can cause skin irritation, burns or blisters. Check your skin regularly.
- Remove the product immediately if you experience skin irritation or if the patch becomes uncomfortably hot.

PRECAUTIONS

- If the contents contact skin or eyes, rinse immediately with water and seek medical advice.
- Do not apply to the same area for more than 12 hours in a 24-hour period.
- If symptoms persist or worsen, consult a doctor.

- Do not exert pressure on the patch for long periods of time (e.g. through the belt or waistband).
- Do not use after drinking alcohol.
- Dispose of in household waste after use.
- Do not cut, rub, fold or moisten the patch. Do not use if the patch has been damaged.

INTENDED PATIENT GROUP

- Children aged 12 and over.
- Adults who can remove the patch themselves.
- In case of pregnancy, diabetes mellitus, weak circulation, rheumatoid arthritis, nerve damage, skin disease or allergy to patch adhesives, a doctor should be consulted.



Unique product identification



Do not reuse



Don't use it when packaging is damaged



Temperature restriction



Expiry date



Medical device



Keep away from sunlight



Keep away from water



Batch designation



CE mark with identification number of the appointed body



Read leaflet



CT & T Medical GmbH
Pfarrmatte 6, CH-8807
Freienbach, Switzerland
www.chinatrading.ch



Sallmann Medical GmbH
Alois-Stockinger-Straße 16
5020 Salzburg, Austria
www.sallmann-medical.eu